	Do not use if package is damaged Niet gebruiken indien de verpakking beschadigd is Ne pas utiliser si l'emballage est endommagé Nicht verwenden, wenn die Verpackung beschädigt ist Non utilizzare se la confezione è danneggiata No lo utilice si el envase está dañado Não usar se a embalagem estiver danificada لا تستعمل التجهيز إذا كانت العبوة تالفة		Do not re-use Niet herbruiken Ne pas réutiliser Nicht wiederverwenden Non riutilizzare No lo reutilice Não reutilizar لا تعد الاستعمال		Do not resterilize Niet opnieuw steriliseren Ne pas restériliser Nicht resterilisieren Non risterilizzare No reesterilizar Não reesterilizar لا تقم بإعادة التعقيم		Keep dry Droog houden Conserver au sec Trocken halten Conservare in un luogo asciutto Manténgase en un lugar seco Mantenha seco أبقه جافاً
--	--	---	---	--	---	---	---

	Manufacturer Producent Fabricant Hersteller Produttore Fabricante المشركة المُصنَّعة.		Manufacturing date Productiedatum Date de fabrication Herstellungsdatum Data di produzione Fecha de fabricación Data de fabricação تاريخ التصنيع		Catalog number Catalogusnummer Référence catalogue Katalognummer Numero di catalogo Número de catálogo Número de catálogo رقم الكتالوج		Lot number Lotnummer Numéro de lot Chargenbezeichnung Numero di lotto Número de lote Número do lote رقم الدفعة		Expiry date Vervaldatum Date de péremption Verfalldatum Data di scadenza Fecha de caducidad Data de expiração تاريخ انتهاء الصلاحية
--	---	---	---	---	---	--	---	---	--

	Fragile, handle with care Breekbaar, met zorg behandelen Fragile, manipuler avec précaution Zerbrechlich – mit Vorsicht handhaben! Materiale fragile, maneggiare con cura Fragil, manejar con cuidado Material frágil, manuseie com cuidado قابل للكسر، تعامل معه بعناية		Keep away from sunlight Weg houden van zonlicht Conserver à l'abri de la lumière du soleil Vor Sonnenlicht schützen Tenere lontano dalla luce del sole Manténgase alejado de la luz solar Mantenha protegido da luz solar يحفظ بعيداً عن ضوء الشمس		Temperature limit Temperatuurlimiet Limite de température Temperaturlimit Limite di temperatura Limite de temperatura Limite de temperatura حد درجة الحرارة
--	---	---	---	---	--

	Unique Device Identifier Unieke identificatiecode voor medische hulpmiddelen Identifiant unique du dispositif Einmalige Produktkennung Identificazione unica del dispositivo Identificador único del producto Identificação única de dispositivo رمز تعريف الجهاز الفريد		CE Marking CE-markering Marquage CE CE-Kennzeichnung Marcatura CE Marcado CE Marcação CE علامة CE		Medical device Medisch hulpmiddel Dispositif médical Medizinprodukt Dispositivo medico Producto sanitario Dispositivo médico جهاز طبي		Non-sterile Niet-steriel Non stérile Unsteril Non sterile No estéril Não estéril غير معقم
--	---	---	--	--	--	---	--

VAX-ID®
VAX-ID® 2200

REF D00037 - VAX-ID® 2200 32G to be sterilized - UDI-DI 05430003573037
D00036 - VAX-ID® 2200 30G to be sterilized - UDI-DI 05430003573051
D00035 - VAX-ID® 2200 27G to be sterilized - UDI-DI 05430003573044

CE
1282

to be sterilized
steriliseren
à stériliser
sterilisieren
da sterilizzare
para ser esterilizado
a ser esterilizado
يجب أن يعقم

English (EN)	Processing instructions	3
	User instructions	4
Nederlands (NL)	Verwerkingsinstructies	5
	Gebruiksaanwijzing	6
Français (FR)	Instructions de traitement	7
	Instructions pour l'utilisateur	8
Deutsch (DE)	Verarbeitungshinweise	9
	Gebrauchsanweisung	10
Italiano (IT)	Istruzioni per la lavorazione	11
	Istruzioni per l'uso	12
Español (ES)	Instrucciones de procesamiento	13
	Instrucciones para el usuario	14
Português (PT)	Instruções de processamento	15
	Instruções de uso	16
عربي (AR)	تعليمات المعالجة	17
	تعليمات للمستخدم	18

Translation
Vertaling
Traduction
Übersetzung
Traduzione
Traducción
Tradução
ترجمة

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
LANDOOR SRL
Via Copernico 38 Milan Italy

IDEVAX BV

Bijkhoevelaan 32c, 2110 Wijnegem, BELGIUM
info@idevax.com www.idevax.com

Cycle phase	
Vacuum phase	Vacuum set point: -80 (+/-5) kPa Temperature raise between pulses: 105 - 110°C Number of pulses: 3
Sterilization phase	High pressure set point: 110 (+5) kPa Temperature: 121°C Duration: 20 min
Drying phase	Low vacuum set point: -40 (+/-5) kPa Temperature: 65 (+/-5)°C Duration: 10 min or more

TABLE 1

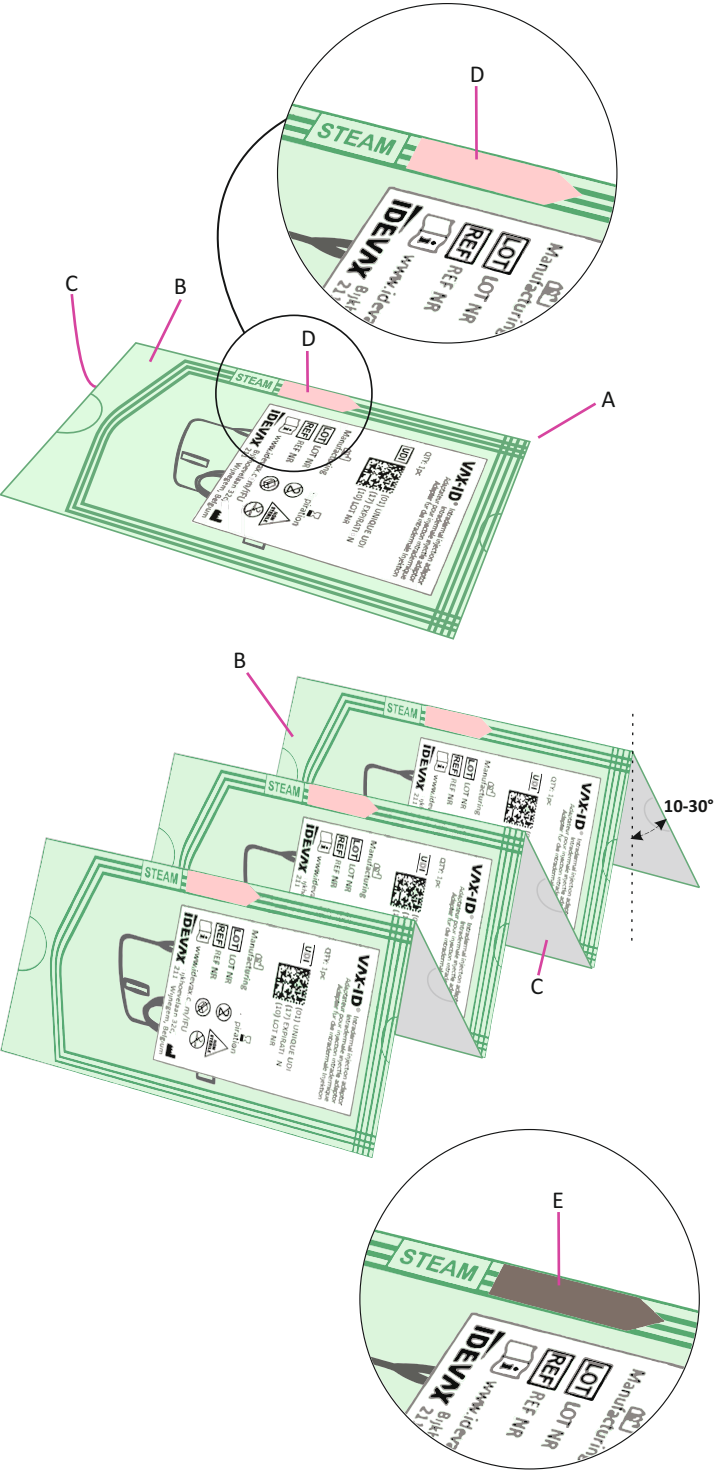EN: Sterilization instructions

DESCRIPTION & INDICATIONS

This instruction provides the Steam (Moist-heat) sterilization process of the VAX-ID® 2200 device in accordance with a validated cycle and EN ISO 17665-1:2024 prior use of the device. This instruction is applicable only for the VAX-ID® variants mentioned in this document. This instruction was generated in compliance with ISO 17664-1:2021

EQUIPMENT AND MATERIALS

- 1. Autoclave
- 2. Autoclave trays

WARNINGS AND PRECAUTIONS

- 1. NEVER open the pouch during any step of this instruction. The pouch is never opened for the sterilization process.
- 2. The device MUST be sterile for use.
- 3. Local or national specifications should be followed where steam sterilization requirements are stricter or more conservative than those listed in table 1.
- 4. The facility doing the sterilization is responsible for in-house procedures for the calibration and maintenance of the autoclave and its work facilities and appropriate personal training.
- 5. VAX-ID® is a single-use device and therefore shall not be sterilized a second time nor reused.

PROCEDURE

The VAX-ID® 2200 device comes in a sterile barrier (A) for the sterilization process. Therefore, the pouch shall not be opened prior to sterilization or use. All pouches have a transparent plastic foil front with label (B) and a paper back (C)

STEP 1: Tray loading

The pouches should be positioned in a way that the paper back is free for moist to enter. Ideally, the pouches are therefore stacked vertically with an angle of 10-30° paper side (C) downwards, in an alternating paper-paper (C) and foil-foil orientation (B). Tray should not be overpacked to allow moist in and to facilitate the drying process at the end of the procedure.

STEP 2: Sterilisation protocol



The steam sterilisation must be performed in a well calibrated and maintained autoclave by trained personnel as per customer’s internal procedures. The sterilisation cycle recommended is described in table 1 and has to be validated by the user, for recommended minimum sterilization parameters that have been validated by Idevax to provide a 10⁻⁶ sterility assurance level (SAL).

STEP 3: Quality control

Check the steam sterilization colour indicator on the pouch which should change from a pink colour (D) before sterilization to a brown colour (E) after sterilization. If the sterilization indicator did not change, do not use the device, and contact Idevax. If you find holes, cuts, or other damage in the pouch, do not use that device, and contact Idevax.

STEP 4: Storage and Transportation

To transport the devices to the point of use or other facility just place the devices inside the same cardboard box(es) they arrived in and close the box using tape. For storage, the storage recommendations on the Instructions for use (IFU) should be followed.

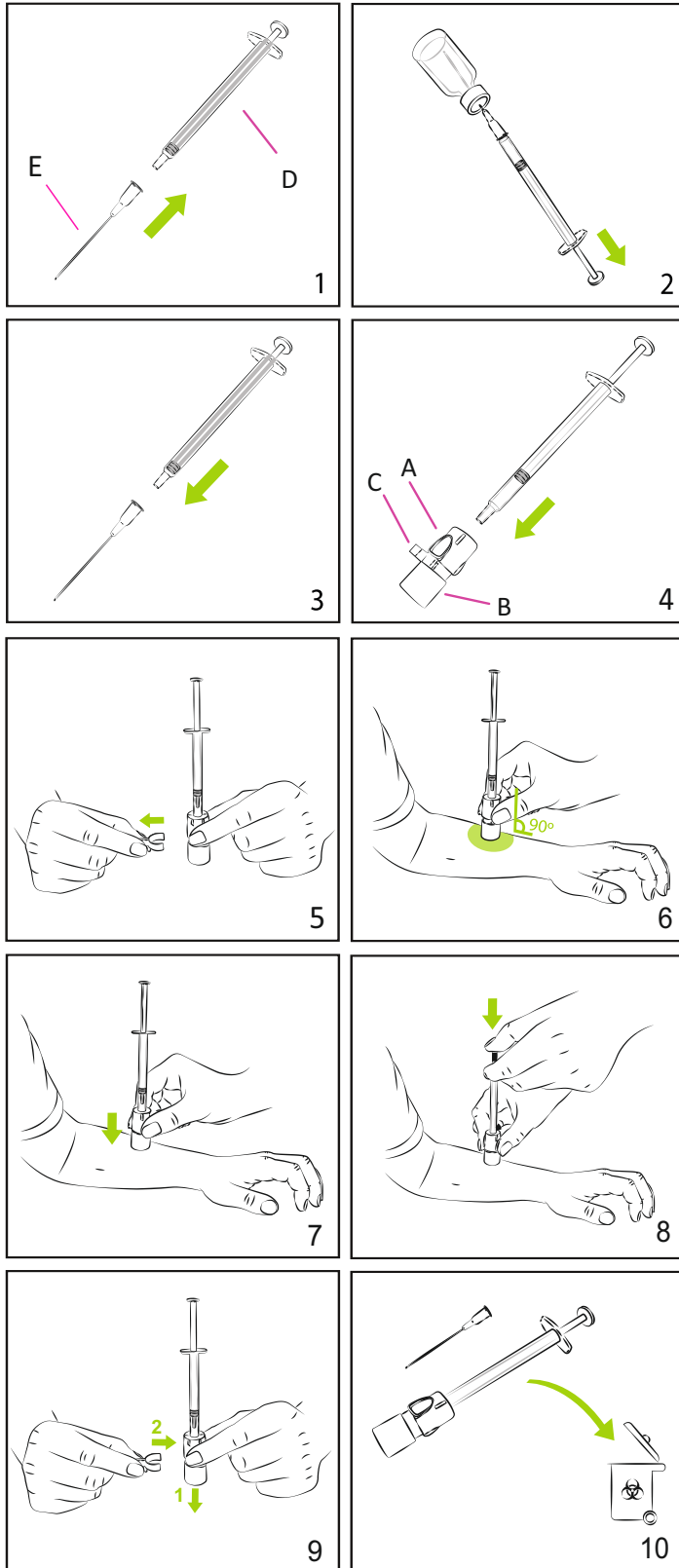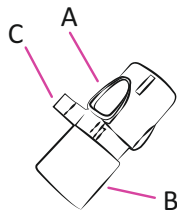EN: User instructions

DESCRIPTION & INDICATIONS

Intradermal injection adaptor: VAX-ID® is a syringe adaptor for easy and reliable intradermal injection of medicinal products. Read all instructions prior to use.

INTENDED USER/ USER QUALIFICATIONS

Healthcare professional (E.g. nurse, doctor, pharmacist)

TARGET POPULATION

All ages and populations.

MATERIALS PROVIDED WITH EACH DEVICE

One adaptor made of:

- A. Housing with injection needle
- B. Foot
- C. Pin

DEVICES TO BE USED IN COMBINATION WITH THE MEDICAL DEVICE

- D. 1.0 mL Luer slip fit syringe (Recommended: Low Dead Space)
- E. Draw needle

WARNINGS AND PRECAUTIONS

Do not use if package is damaged or open.

Use only after sterilization.

Single Use Only. Re-use may lead to contamination by foreign material/organisms thus infection and/or illness.

Dispose of used device in a sharp's container in accordance with applicable national regulations.

Any serious incidence (as defined per the European Medical Device Regulation 2017/745) should be reported to the manufacturer and the competent Authority.

ADVERSE EVENTS

Intradermal injections may be associated with transient local skin reactions including local pain, edema, erythema, itching, hardening, peeling, discoloration at the injection site, most likely related to the injected substance, and not the device.

Adverse events should be reported to the manufacturer.

CONTRAINDICATIONS

Do not use on wounds, scars, nor damaged skin.

STORAGE

The device should be stored at temperatures ranging from 10-30°C, in a dry environment, with no direct sunlight, and without large temperature fluctuations.

INSTRUCTIONS FOR USE

1. Attach the draw needle (E) to the syringe (D).
2. Draw the volume to be injected from the vial.
3. Remove the draw needle from the syringe.
4. Put the syringe into VAX-ID® adaptor via the opening on the end of the housing (A). Make sure that the syringe tip is properly fitted into the Luer fitting of the embedded injection needle.
5. Activate VAX-ID® by removing the pin (C).
6. Place the foot (B) on the skin, keeping the adaptor in a perpendicular position.
7. Push the housing (A) down to meet the foot (B): the needle will penetrate the skin.
8. Retain slight pressure with two fingers on the housing (A) and use another finger to slowly (3-5 seconds) inject the fluid by pushing the syringe plunger, expect resistance while injecting, and wait few (2-3) seconds before removing the VAX-ID®.
9. Remove VAX-ID® from the skin and deactivate the device by pulling the foot (B) away from the housing (A) and placing the pin (C) back in its original position.
10. Dispose VAX-ID® and the draw needle in a sharps container.

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5. Aktivieren Sie VAX-ID®, indem Sie den Stift © entfernen.
VAX-ID is a registered trademark in the name of Idevax BV / VAX-ID is een geregistreerd handelsmerk van Idevax BV