|                                                                                  |                                                                                                                                                                                                                                                                                                                                                                                  |                                                                                   |                                                                                                                                                             |                                                                                    |                                                                                                                                                                                         |                                                                                     |                                                                                                                                                                 |
|----------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | Do not use if package is damaged<br>Niet gebruiken indien de verpakking beschadigd is<br>Ne pas utiliser si l'emballage est endommagé<br>Nicht verwenden, wenn die Verpackung beschädigt ist<br>Non utilizzare se la confezione è danneggiata<br>No lo utilice si el envase está dañado<br>Não usar se a embalagem estiver danificada<br>لا تستعمل التجهيز إذا كانت العبوة تالفة |  | Do not re-use<br>Niet herbruiken<br>Ne pas réutiliser<br>Nicht wiederverwenden<br>Non riutilizzare<br>No lo reutilice<br>Não reutilizar<br>لا تعد الاستعمال |  | Do not resterilize<br>Niet opnieuw steriliseren<br>Ne pas restériliser<br>Nicht resterilisieren<br>Non risterilizzare<br>No reesterilizar<br>Não reesterilizar<br>لا تقم بإعادة التعقيم |  | Keep dry<br>Droog houden<br>Conserver au sec<br>Trocken halten<br>Conservare in un luogo asciutto<br>Manténgase en un lugar seco<br>Mantenha seco<br>أبقه جافاً |
|----------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------|

|                                                                                  |                                                                                                         |                                                                                   |                                                                                                                                                                       |                                                                                   |                                                                                                                                                             |                                                                                    |                                                                                                                                     |                                                                                     |                                                                                                                                                          |
|----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | Manufacturer<br>Producent<br>Fabricant<br>Hersteller<br>Produttore<br>Fabricante<br>المشركة المُصنَّعة. |  | Manufacturing date<br>Productiedatum<br>Date de fabrication<br>Herstellungsdatum<br>Data di produzione<br>Fecha de fabricación<br>Data de fabricação<br>تاريخ التصنيع |  | Catalog number<br>Catalogusnummer<br>Référence catalogue<br>Katalognummer<br>Numero di catalogo<br>Número de catálogo<br>Número de catálogo<br>رقم الكتالوج |  | Lot number<br>Lotnummer<br>Numéro de lot<br>Chargenbezeichnung<br>Numero di lotto<br>Número de lote<br>Número do lote<br>رقم الدفعة |  | Expiry date<br>Vervaldatum<br>Date de péremption<br>Verfalldatum<br>Data di scadenza<br>Fecha de caducidad<br>Data de expiração<br>تاريخ انتهاء الصلاحية |
|----------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------|

|                                                                                  |                                                                                                                                                                                                                                                                                               |                                                                                   |                                                                                                                                                                                                                                                                         |                                                                                     |                                                                                                                                                                                  |
|----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | Fragile, handle with care<br>Breekbaar, met zorg behandelen<br>Fragile, manipuler avec précaution<br>Zerbrechlich – mit Vorsicht handhaben!<br>Materiale fragile, maneggiare con cura<br>Frágil, manejar con cuidado<br>Material frágil, manuseie com cuidado<br>قابل للكسر، تعامل معه بعناية |  | Keep away from sunlight<br>Weg houden van zonlicht<br>Conserver à l'abri de la lumière du soleil<br>Vor Sonnenlicht schützen<br>Tenere lontano dalla luce del sole<br>Manténgase alejado de la luz solar<br>Mantenha protegido da luz solar<br>يحفظ بعيداً عن ضوء الشمس |  | Temperature limit<br>Temperatuurlimiet<br>Limite de température<br>Temperaturlimit<br>Limite di temperatura<br>Límite de temperatura<br>Limite de temperatura<br>حد درجة الحرارة |
|----------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|

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LANDOOR SRL  
Via Copernico 38 Milan Italy

**IDEVAX BV**

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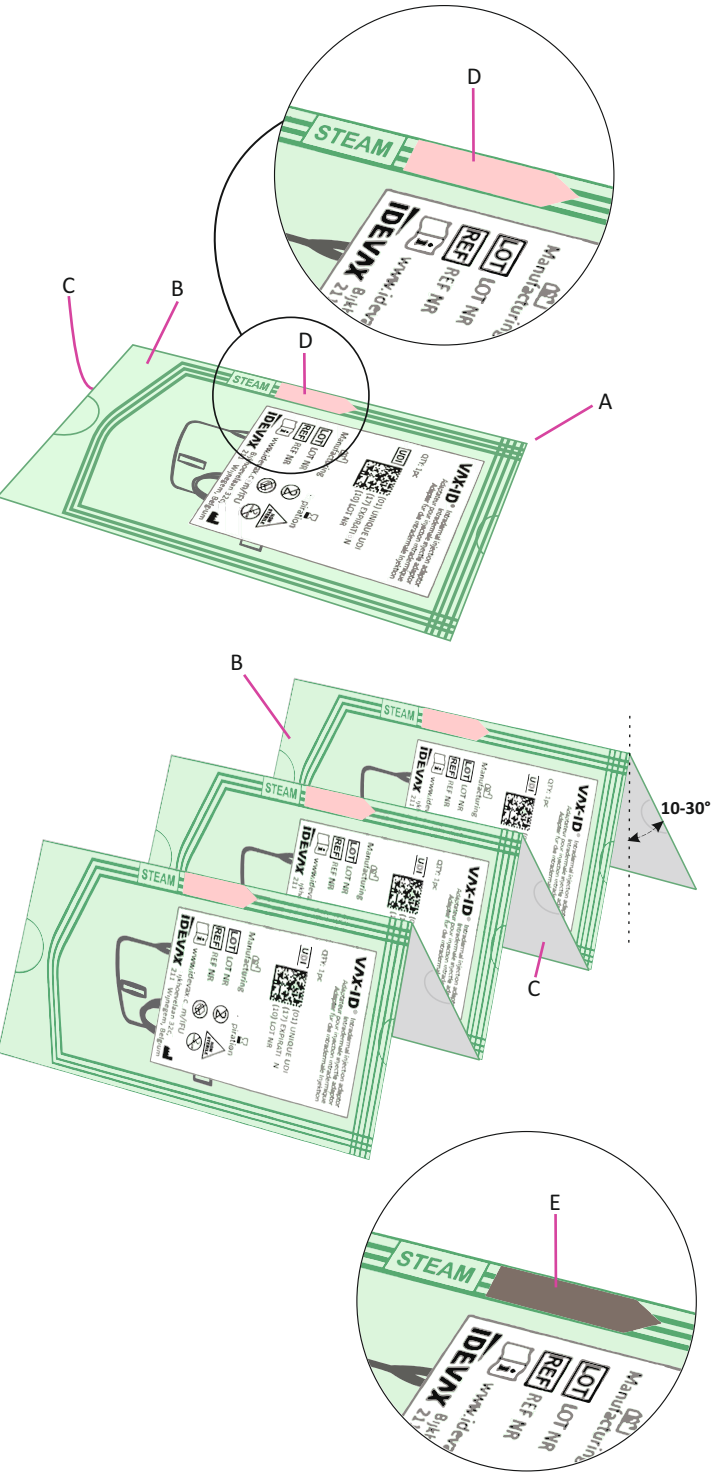
[info@idevax.com](mailto:info@idevax.com) [www.idevax.com](http://www.idevax.com)

| Cycle phase         |                                                                                                                               |
|---------------------|-------------------------------------------------------------------------------------------------------------------------------|
| Vacuum phase        | Vacuum set point: <b>-80 (+/-5) kPa</b><br>Temperature raise between pulses: <b>105 - 110°C</b><br>Number of pulses: <b>3</b> |
| Sterilization phase | High pressure set point: <b>110 (+5) kPa</b><br>Temperature: <b>121°C</b><br>Duration: <b>20 min</b>                          |
| Drying phase        | Low vacuum set point: <b>-40 (+/-5) kPa</b><br>Temperature: <b>65 (+/-5)°C</b><br>Duration: <b>10 min or more</b>             |

TABLE 1

EN: Sterilization instructions

DESCRIPTION & INDICATIONS

This instruction provides the Steam (Moist-heat) sterilization process of the VAX-ID® 2200 device in accordance with a validated cycle and EN ISO 17665-1:2024 prior use of the device. This instruction is applicable only for the VAX-ID® variants mentioned in this document. This instruction was generated in compliance with ISO 17664-1:2021

EQUIPMENT AND MATERIALS

- 1. Autoclave
- 2. Autoclave trays

WARNINGS AND PRECAUTIONS

- 1. NEVER open the pouch during any step of this instruction. The pouch is never opened for the sterilization process.
- 2. The device MUST be sterile for use.
- 3. Local or national specifications should be followed where steam sterilization requirements are stricter or more conservative than those listed in table 1.
- 4. The facility doing the sterilization is responsible for in-house procedures for the calibration and maintenance of the autoclave and its work facilities and appropriate personal training.
- 5. VAX-ID® is a single-use device and therefore shall not be sterilized a second time nor reused.

PROCEDURE

The VAX-ID® 2200 device comes in a sterile barrier (A) for the sterilization process. Therefore, the pouch shall not be opened prior to sterilization or use. All pouches have a transparent plastic foil front with label (B) and a paper back (C)

STEP 1: Tray loading

The pouches should be positioned in a way that the paper back is free for moist to enter. Ideally, the pouches are therefore stacked vertically with an angle of 10-30° paper side (C) downwards, in an alternating paper-paper (C) and foil-foil orientation (B). Tray should not be overpacked to allow moist in and to facilitate the drying process at the end of the procedure.

STEP 2: Sterilisation protocol

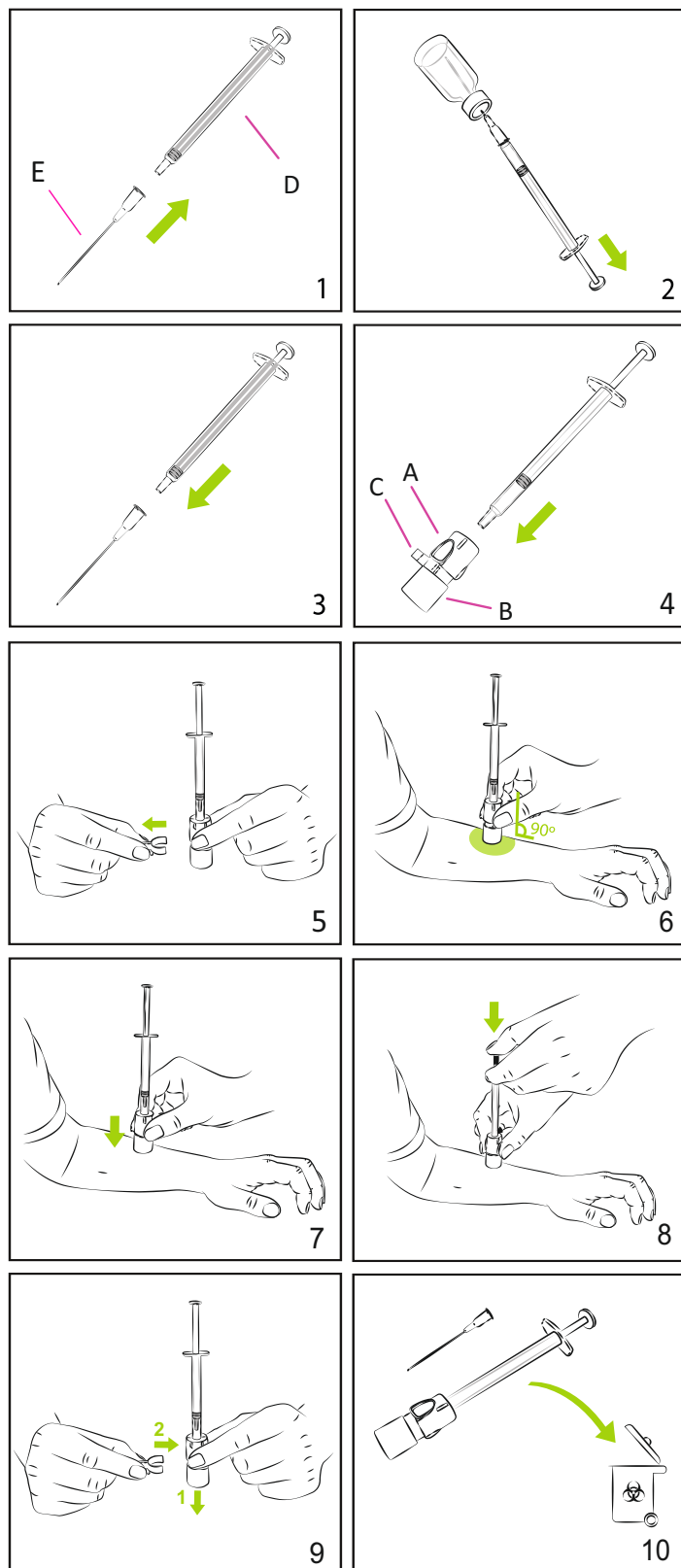
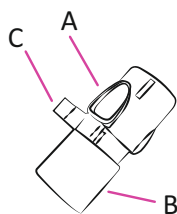
The steam sterilisation must be performed in a well calibrated and maintained autoclave by trained personnel as per customer’s internal procedures. The sterilisation cycle recommended is described in table 1 and has to be validated by the user, for recommended minimum sterilization parameters that have been validated by Idevax to provide a 10<sup>-6</sup> sterility assurance level (SAL).

STEP 3: Quality control

Check the steam sterilization colour indicator on the pouch which should change from a pink colour (D) before sterilization to a brown colour (E) after sterilization. If the sterilization indicator did not change, do not use the device, and contact Idevax. If you find holes, cuts, or other damage in the pouch, do not use that device, and contact Idevax.

STEP 4: Storage and Transportation

To transport the devices to the point of use or other facility just place the devices inside the same cardboard box(es) they arrived in and close the box using tape. For storage, the storage recommendations on the Instructions for use (IFU) should be followed.

## EN: User instructions

### DESCRIPTION & INDICATIONS

Intradermal injection adaptor: VAX-ID® is a syringe adaptor for easy and reliable intradermal injection of medicinal products.

Read all instructions prior to use.

### INTENDED USER/ USER QUALIFICATIONS

Healthcare professional (E.g. nurse, doctor, pharmacist)

### TARGET POPULATION

All ages and populations.

### MATERIALS PROVIDED WITH EACH DEVICE

One adaptor made of:

A. Housing with injection needle

B. Foot

C. Pin

### DEVICES TO BE USED IN COMBINATION WITH THE MEDICAL DEVICE

D. 1.0 mL Luer slip fit syringe (Recommended: Low Dead Space)

E. Draw needle

### WARNINGS AND PRECAUTIONS

Do not use if package is damaged or open.

Use only after sterilization.

Single Use Only. Re-use may lead to contamination by foreign material/organisms thus infection and/or illness.

Dispose of used device in a sharps container in accordance with applicable national regulations.

Any serious incidence (as defined per the European Medical Device Regulation 2017/745) should be reported to the manufacturer and the competent Authority.

### ADVERSE EVENTS

Intradermal injections may be associated with transient local skin reactions including local pain, edema, erythema, itching, hardening, peeling, discoloration at the injection site, most likely related to the injected substance, and not the device.

Adverse events should be reported to the manufacturer.

### CONTRAINDICATIONS

Do not use on wounds, scars, nor damaged skin.

### STORAGE

The device should be stored at temperatures ranging from 10-30°C, in a dry environment, with no direct sunlight, and without large temperature fluctuations.

### INSTRUCTIONS FOR USE

1. Attach the draw needle (E) to the syringe (D).
2. Draw the volume to be injected from the vial.
3. Remove the draw needle from the syringe.
4. Put the syringe into VAX-ID® adaptor via the opening on the end of the housing (A). Make sure that the syringe tip is properly fitted into the Luer fitting of the embedded injection needle.
5. Activate VAX-ID® by removing the pin (C).
6. Place the foot (B) on the skin, keeping the adaptor in a perpendicular position.
7. Push the housing (A) down to meet the foot (B): the needle will penetrate the skin.
8. Retain slight pressure with two fingers on the housing (A) and use another finger to slowly (3-5 seconds) inject the fluid by pushing the syringe plunger, expect resistance while injecting, and wait few (2-3) seconds before removing the VAX-ID®.
9. Remove VAX-ID® from the skin and deactivate the device by pulling the foot (B) away from the housing (A) and placing the pin (C) back in its original position.
10. Dispose VAX-ID® and the draw needle in a sharps container.

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5. Aktivieren Sie VAX-ID®, indem Sie den Stift © entfernen.